

**KARMA**  
EXPORTS

# PROTECTIVE MASK GB 2626-2006 KN95



MORE  
PROTECTIVE  
  
MORE  
COMFORTABLE



web : [www.karmaexports.in](http://www.karmaexports.in)



# PROTECTIVE MASK GB 2626-2006 KN95

## MORE PROTECTIVE

### FLEXIBLE NOSEBAND

Adjustable clip makes the Mask fit your face better.

### HIGH STRENGTH ELASTICITY EARLOOPS

Normal use will not easily break, but do not pull very hard.

Adjustable nose clip

Breathable and Comfortable

Enhance Protection


Four-Layer Design

① PP non-woven

② Meltblown fabric

③ Meltblown fabric


④ PP non-woven

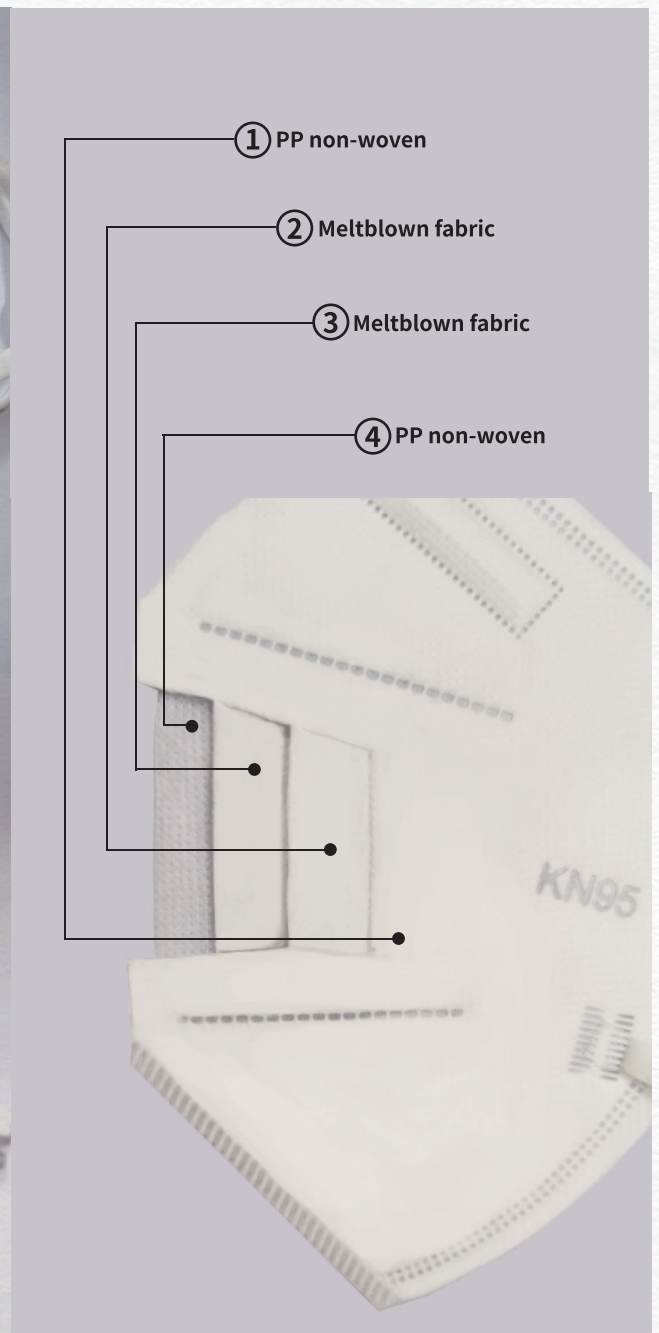



## 50 Pieces/Pack with LOGO Individual package

**Size:** 160\*125mm  
**Material:** OPP+CPP

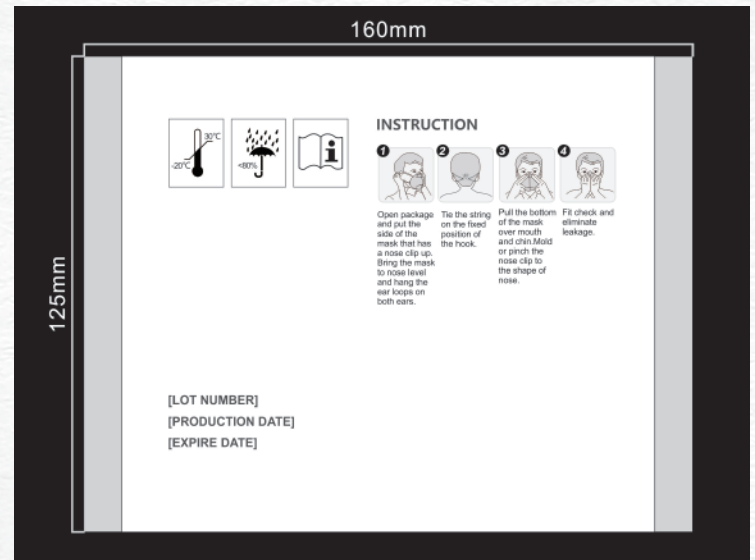

## Color Box

**Size:** 205\*155\*128mm  
**Material:** 400g Film Lamination



## Master Box

**Size:** 650\*425\*540mm





# QUALIFIED CERTIFICATE

**NAME:**

Protective Mask

**MODEL:**

Ear hanging type (KF-A F04)

**STANDARD:**

GB 2626-2006

**MAIN COMPONENTS:**

PP nonwoven fabric, PP melt-blown fabric

**PRODUCTION DATE:**

Printed on the package

**LOT NUMBER:**

Printed on the package

**VALIDITY PERIOD:**

2 Years





# INSTRUCTION MANUAL

---

**[Name]** Protective Mask

**[Model]** Ear hanging type (KF-A F04)

**[Applicable standard]** GB 2626-2006

**[Specification]** 23cm×12cm

**[Structure and composition]** This product consists of PP nonwoven fabric, PP melt-blown fabric, nose clip and ear loops.

**[Intended use]**

The protective mask is intended to be worn to protect personnel from the transfer of particulate material. This is a single use, non-sterile product.

**[Instructions for use]**

1. Open package and put the side of the mask that has a nose clip up. Bring the mask to nose level and hang the ear loops on both ears;
2. Tie the string on the fixed position of the hook;
3. Pull the bottom of the mask over mouth and chin. Mold or pinch the nose clip to the shape of nose;
4. Fit check and eliminate leakage.

**[Precautions]**

1. This product is a non-sterile type mask.
2. This product's validity period is 2 years. Please use within validity period.
3. The product is disposable and cannot be reused or recycled. After use, it shall be immediately destroyed or thrown into a special treatment box.
4. Please check the integrity of the package before use. Do not use the product if the package is damaged.
5. Do not use the product if you are allergic to it.
6. Incorrect wearing of this product may cause breathing difficulties.
7. Children under three years old are not recommended to use the product because of their low vital capacity.
8. The mask should be kept as flat as possible after use, and the nose clip above the mask should not be folded frequently to prolong the service time.
9. Used masks are considered highly contaminated, it is essential that:
  - a) The body of the mask is not touched by the fingers/hands of the wearer;
  - b) Hands are disinfected (full hand disinfection) after mask removal;
  - c) A mask is worn covering the nose and mouth of the wearer, at no time a mask is hanging around the neck of the wearer.

**[Validity period]** 2 years

**[Storage and transportation conditions]** The packaged mask shall be stored in a room with normal temperature and relative humidity of no more than 80%, no corrosive gas, good ventilation, and fire, rat and insect prevention facilities. During transportation, keep away from moisture, light and heat.

**[Package]** 50pcs/bag.












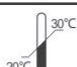
**[Production date]** See package

**[Lot number]** See package



# INSTRUCTION MANUAL

## [GRAPHS AND SYMBOLS DESCRIPTION]

1		See information supplied by the manufacturer
2		Do not use if the package is damaged
3		Do not re-use
4		Do not dispose of this product together with domestic household wastes
5		Keep dry
6		This side up
7		Fragile, handle with care
8		Keep away from sunlight
9		Use-by date
10		Lot number
11		Maximum number of stack levels
12		Temperature range of storage conditions



# FACTORY ENVIRONMENT





.....

# KARMA EXPORTS

50, Maddox Street, Near Veterinary College,  
Vepery, Chennai - 600 007, South India.

e mail : [info@karmaexports.in](mailto:info@karmaexports.in)

Tel : +91 44 2661 0819

+91 97908 40995, +91 7708370525, +91 86759 99222

web : [www.karmaexports.in](http://www.karmaexports.in)

.....