

DISPOSABLE PROTECTIVE MASK



MORE
PROTECTIVE

MORE
COMFORTABLE

DISPOSABLE PROTECTIVE MASK

MORE PROTECTIVE

FLEXIBLE NOSEBAND

Adjustable clip makes the Mask fit your face better.



HIGH STRENGTH ELASTICITY EARLOOPS


Normal use will not easily break, but do not pull very hard.

KINGFA

① PP Nonwoven

② Meltblown Fabric

③ PP Nonwoven

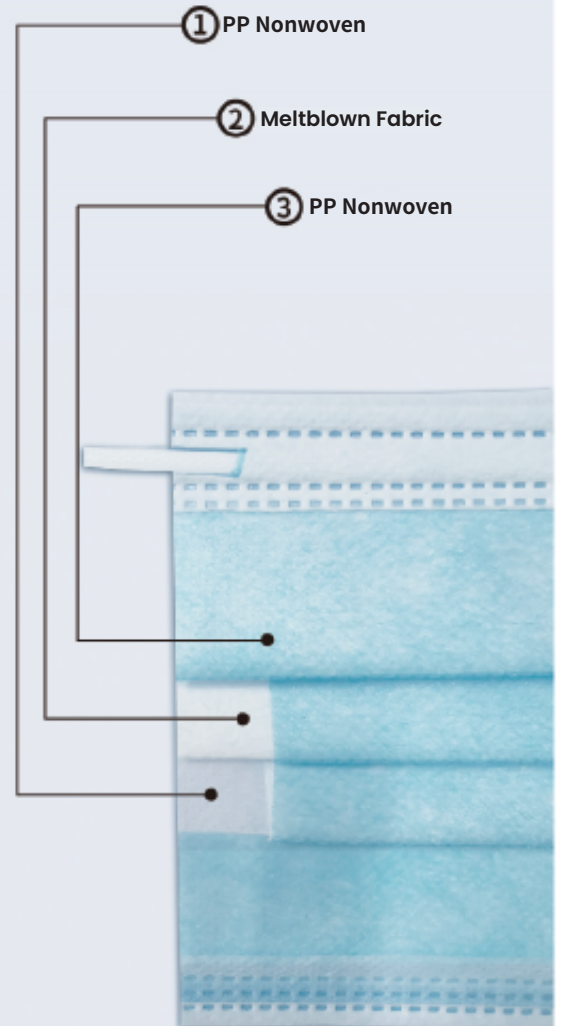


Three-Layer Design

Adjustable nose clip

Breathable and Comfortable


Enhance Protection

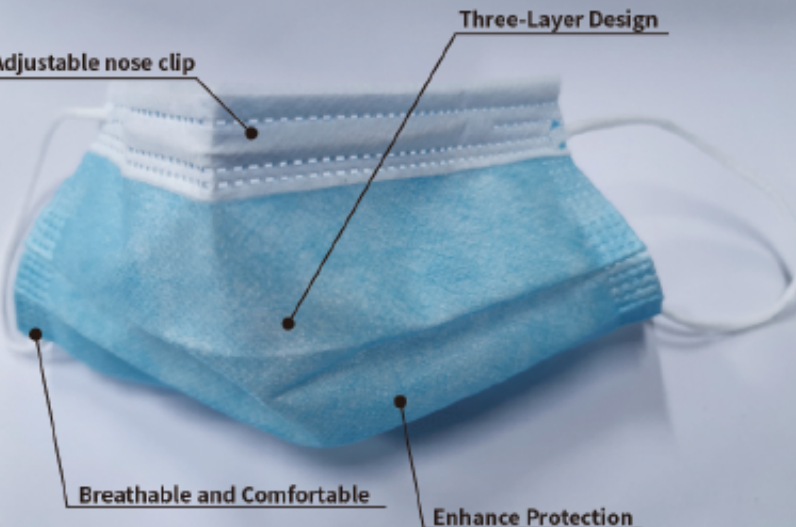

50 Pieces/Pack with LOGO Heat-sealed package

Size: 360*120mm

Material: LDPE

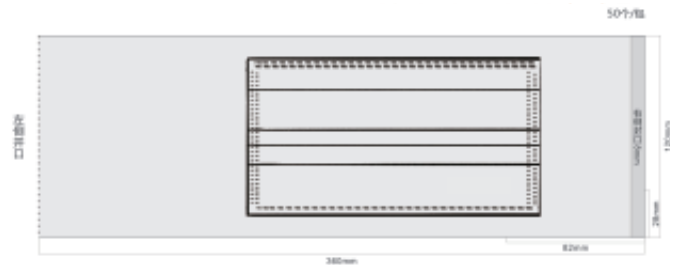

Color Box (50 Pieces/Box)

Size: 210*100*100mm

Material: 400g Film Lamination



Master Box

Size: 530*440*530mm



QUALIFIED CERTIFICATE

NAME:

Disposable Protective Mask

MODEL:

Adult type (KF-B P01)

TRADEMARK:

KF-B P01(TB)

STANDARD:

T/GDBX 025-2020

MAIN COMPONENTS:

PP nonwoven fabric, PP melt-blown fabric

PRODUCTION DATE:

Printed on the package

LOT NUMBER:

Printed on the package

VALIDITY PERIOD:

2 Years



INSTRUCTION MANUAL

[Model] Adult type (KF-B P01)

[Trademark] KF-B P01(TB)

[Specification] Flat type 17.5cm×9.5cm

[Structure and composition] This product consists of PP nonwoven fabric, PPMelt-blown fabric, nose clip and ear loops.

[Intended use]

The disposable protective mask is intended to be worn to protect personnel from the transfer of particulate material. This is a single use, non-sterile product.

[Instructions for use]

1. Open the sealed package and take out the mask.
2. Attach the mask to the face and nose horizontally, with the nose clip pointing upwards.
3. Place the fingertips of both hands on the nose clip, from the middle position, with the fingers to press inward, and gradually move to the sides to make the nose clip of the mask fit the shape of the nose.
4. Pull the ear loops and wear the mask on the ear;
5. Spread the pleated creases and extend the mask.

[Precautions]

1. This product is a non-sterile type mask.
2. This product's validity period is 2 years. Please use within validity period.
3. The product is disposable and cannot be reused or recycled. After use, it shall be immediately destroyed or thrown into a special treatment box.
4. Please check the integrity of the package before use. Do not use the product if the package is damaged.
5. Do not use the product if you are allergic to it.
6. Incorrect wearing of this product may cause breathing difficulties.
7. Children under three years old are not recommended to use the product because of their low vital capacity.
8. The mask should be kept as flat as possible after use, and the nose clip above the mask should not be folded frequently to prolong the service time.
9. Used masks are considered highly contaminated, it is essential that:
 - a) The body of the mask is not touched by the fingers/hands of the wearer;
 - b) Hands are disinfected (full hand disinfection) after mask removal;
 - c) A mask is worn covering the nose and mouth of the wearer, at no time a mask is hanging around the neck of the wearer.

[Validity period] 2 years

[Storage and transportation conditions] The packaged mask shall be stored in a room with normal temperature and relative humidity of no more than 80%, no corrosive gas, good ventilation, and fire, rat and insect prevention facilities. During transportation, keep away from moisture, light and heat.











[Package] 50pcs/bag.

[Production date] See package

[Lot number] See package

INSTRUCTION MANUAL

[GRAPHS AND SYMBOLS DESCRIPTION]

1		Refer to instruction/manual before use
2		Do not use if the package is damaged
3		Do not re-use
4		Do not dispose of this product together with domestic household wastes
5		Keep dry
6		This side up
7		Fragile, handle with care
8		Keep away from sunlight
9		Use-by date
10		Lot number
11		Maximum number of stack levels

FACTORY ENVIRONMENT



.....

KARMA EXPORTS

50, Maddox Street, Near Veterinary College,
Vepery, Chennai - 600 007, South India.

e mail : info@karmaexports.in

Tel : +91 44 2661 0819

+91 97908 40995, +91 7708370525, +91 86759 99222

web : www.karmaexports.in

.....